

Monitor C18 HPLC Column

HPLC Condition:

Column: Monitor C18 HPLC column 2.1x50 mm, 3 μm.

Part number: M3BE-001

Mobile phase:




Time (min)	Water	Acetonitrile	1% Formic Acid
0	85	5	10
7	80	10	10
18	45	45	10

Flow rate: 0.3 mL/min

Injection volume: 10 μL

Detection: UV @ 260 nm



Compounds:

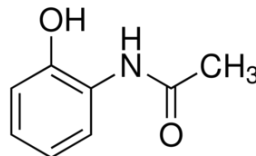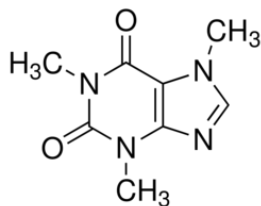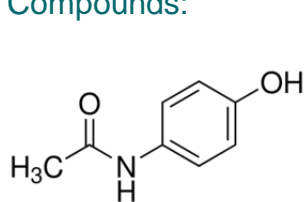

1. Acetaminophen

2. Caffeine

3. 2-Acetamidophenol

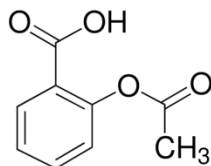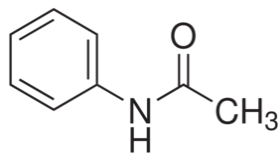

4. Acetanilide

5. Acetylsalicylic Acid

