

Orosil C18 HPLC Column

HPLC Conditions:

Column: Orosil C18 HPLC column 2.1x50 mm, 3 μ m.

Part number: O3BE-902

Mobile phase:


Time (min)	Water	Acetonitrile	1% Formic Acid
0	85	5	10
7	80	10	10
18	45	45	10

Flow rate: 0.3 mL/min


Injection volume: 10 μ L

Detection: UV at 260 nm


Compounds:

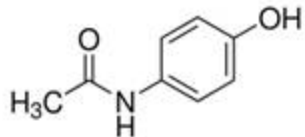


1. Acetaminophen

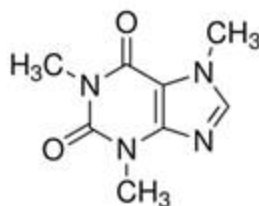


2. Caffeine

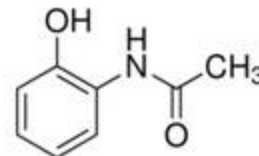

3. 2-Acetamidophenol

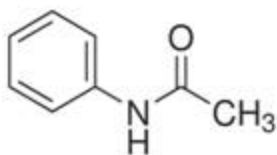

4. Acetanilide

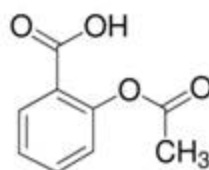

5. Acetylsalicylic Acid

