


Napropamid on Epitomize CSP-1Z, 3 μm , SFC-UHPLC

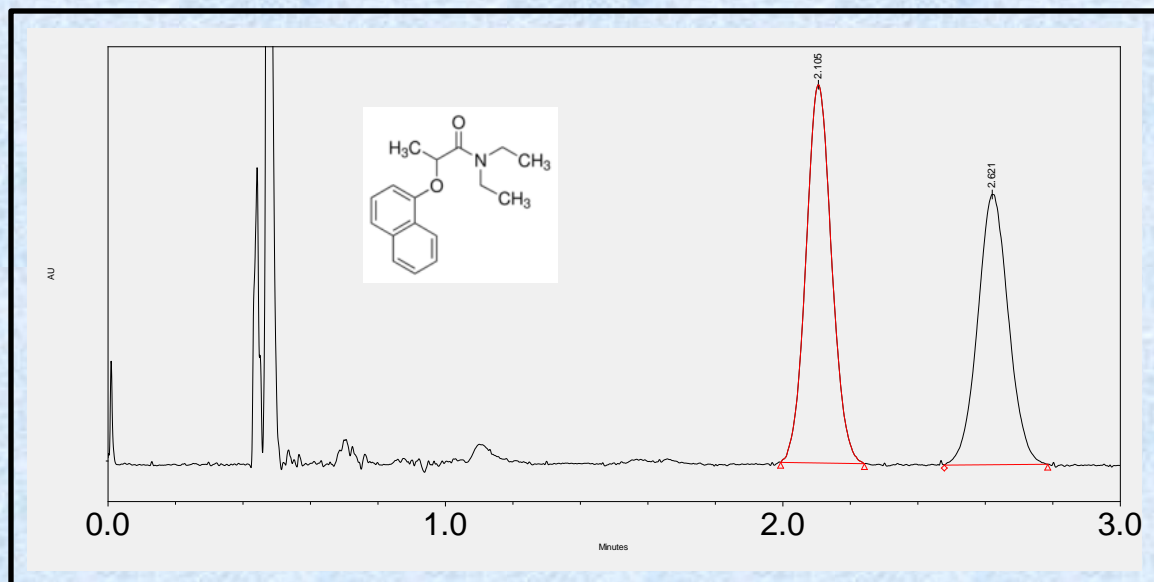

System	Waters Acquity UPC ²
Column	CSP-1Z, 3 micron
Column Size	3.0 x 50 mm
Mobile Phase	90/10 CO ₂ /Modifier
Modifier	50/50 MeOH/IPA
Flow Rate	0.80 mL/min
Temperature	30°C
Detection	UV at 254 nm

Injection	1.0 μL (1 mg/mL)
Rt1 & Rt2 (min)	2.11 & 2.62
k1 & k2	3.40 & 4.48
Separation Value	1.32
Resolution Value	3.39